


Model: SCOTE-36300

Instruction Manual

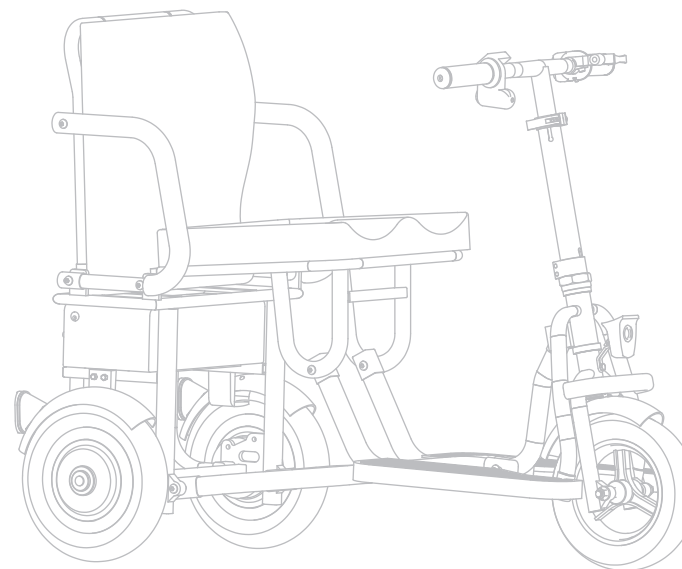

Thank you for your patronage !


Proper use of this product is the prerequisite to ensure safety, convenience and durability

Please read this instruction carefully before using, and be sure to master the matters marked with the mark ★

Catalogue

- 1.Product drawing
- 2.Install
- 3.Opening and folding indication picture
- 4.Accessories description
- 5.Parameters
- 6.Product standard configuration
- 7.Inspection and maintenance before driving
- 8.Battery instruction
- 9.Attentions for safe driving
- 10.Eliminate of frequent malfunction
- 11.Friendly reminder


Product Drawing

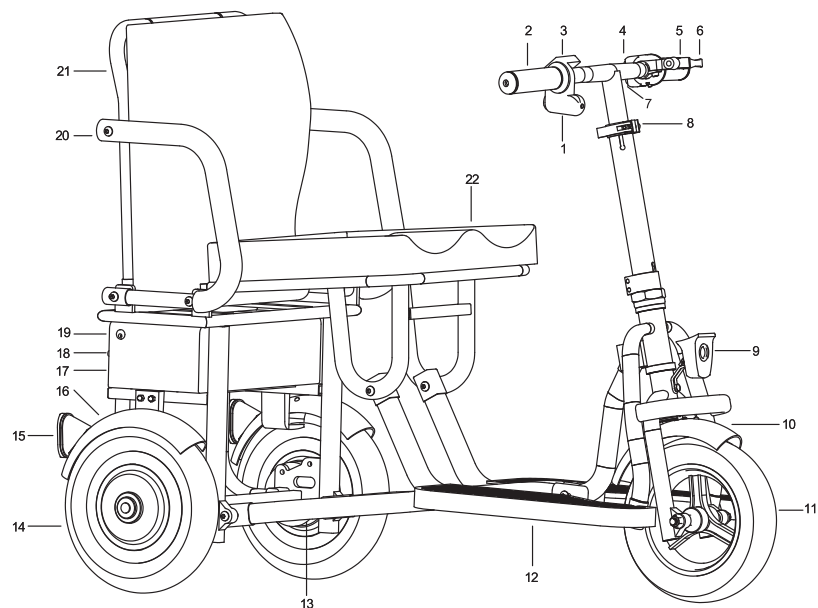


1. Key switch
2. Shifter
3. Electric quantity display
4. Forward and reversed switch
5. Handlebar
6. Brake
7. Horn switch
8. Quick release
9. Lamp
10. Front fender
11. Front wheel

12. Padel plate
13. Brake apparatus
14. Motor
15. Reflector
16. Rear fender
17. Battery
18. Battery box
19. Charge port
20. Armrest
21. Backrest
22. Cushion


Install

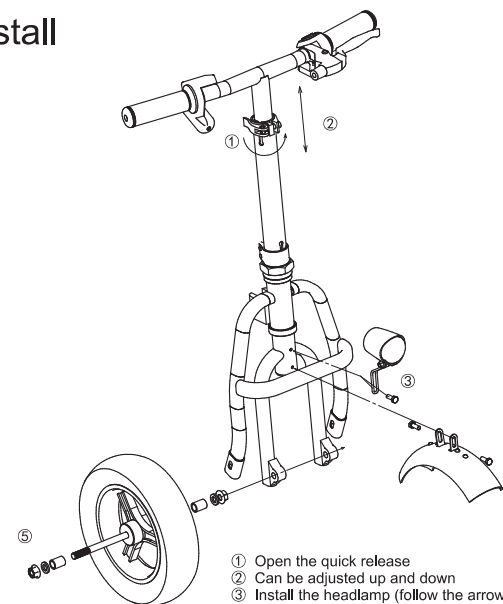

- ① Open the quick release
- ② Can be adjusted up and down
- ③ Install the headlamp (follow the arrow)
- ④ Install the fender (follow the arrow)
- ⑤ Install the wheel (follow the arrow), the wheel should be in the centre

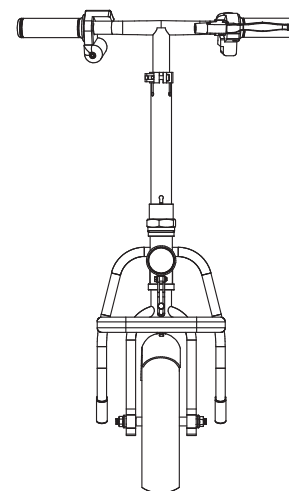


Face side install finish

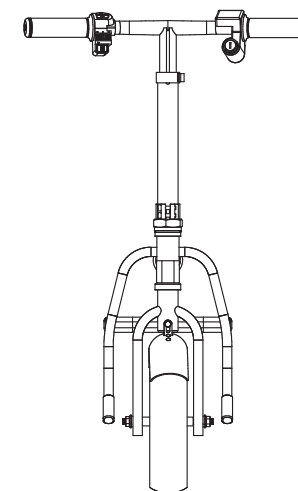


Back side install finish

Opening indication picture

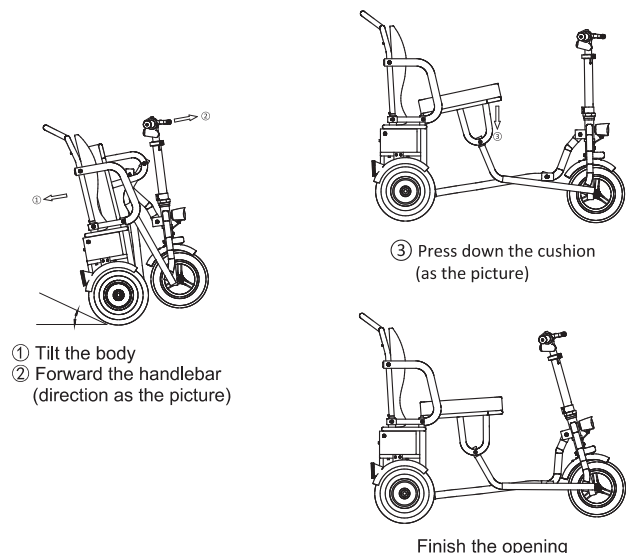


Folding indication picture

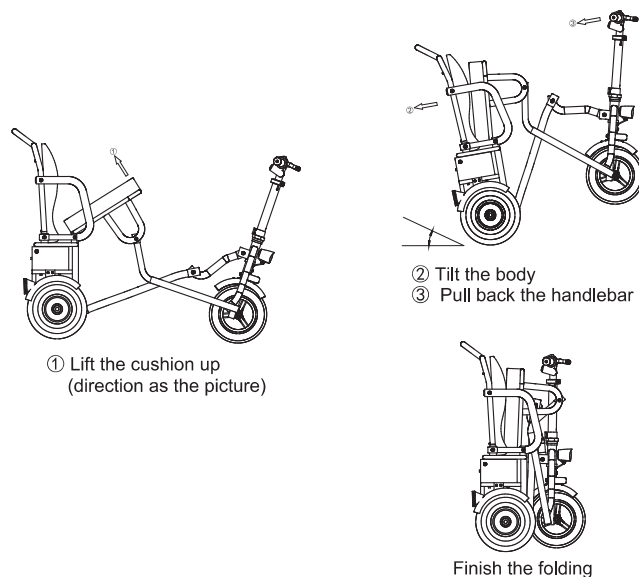

Accessories description

Waterproof plug wire

Waterproof plug, waterproof connector, waterproof performance, can provide safe and reliable electrical, signal and other connections, easy to replace products

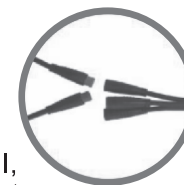

Electric control switch

Electric switch headlights, horns, forward and back switch buttons are designed to be anti-skid, with better symmetry and convenient operation



Red button: forward and reverse shift button

Green button: Lamp switch

Blackout brake handle

Aluminum alloy blackout brake handle, material relative light, not easy to rust, not deformation and brake line coordination can timely brake



Electron turning handle

We focus on the comfortable riding experience. The grip we use is more in line with the manual work of the turn handle, and has good skid resistance and a strong sense of control. Display the riding data, so that the rider has a basis and a bottom in mind



Open the key

The first: close the electric

The second: open the electric

The third: open the electric and lamp



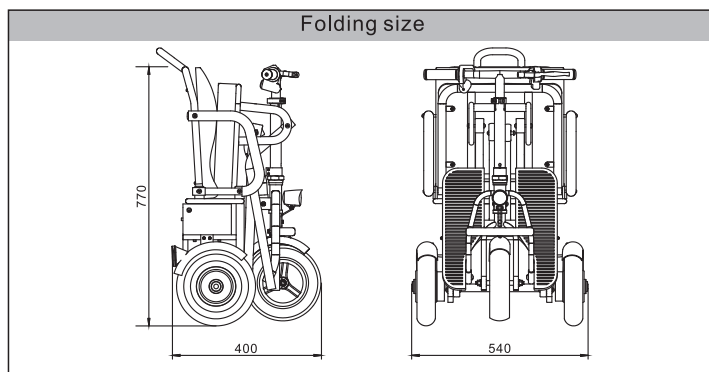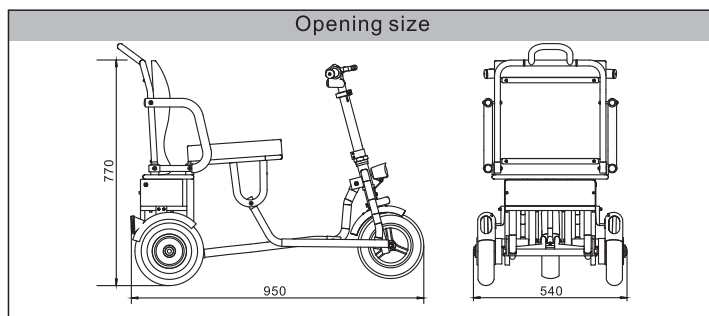

SCOTE-36300 Parameters

Dimensions (mm)				
	Length	Width	Height	
Open	950	540	770	
Fold	400	540	770	
Weight (Kg)		Tyre (Inch)		
G. W	Weight		Front	Back
20.5	≤100		8inch	8inch
Electric specification				
Motor		Lithium battery		
Voltage	Watt	Capacity	Voltage	
36V	300W	8AH	220V	
Charger		Controller		
Voltage	Powder	Current output	Voltage	Current output
36V	220V50HZ	42V-2.0A	36V	15A
specification				
Max. speed	≤15km/h	Back speed	≤8km/h	
Driving	≤30Km	Braking	2500mm	

Product standard configuration

Portable electric folding scooter	Installation tool	Charger	Instructions
1 unit	1unit	1unit	1unit

The above data are for reference only. The technical data of this product series are subject to change without prior notice. In actual use, due to different road conditions, climate and driving speed, data such as charging time and driving range are different.

Inspection and maintenance before driving

In order to ensure that this portable electric folding car can bring convenience to your life, please check and maintain the product according to the following requirements before using.

1. Check whether the battery is full of electricity to ensure that all kinds of cables are properly connected.
2. Check whether the brake system is sensitive and effective.
3. Check whether the tire pressure is normal.
4. Before driving, please correct the front wheel and tighten the handle fastening handle.
5. The handlebar and the moving clearance of each part, should be inspected fastening and lubricated irregularly.

Battery Instruction

★Improper use will shorten the service life of the battery, and may lead to the fever of the battery, Fire, or even explosions.

Please read the following notes carefully:

★Do not separate and decompose the battery. There are protective circuit function parts inside the battery pack.

Disassembly will damage these parts;

★It is strictly prohibited to short circuit the positive and negative terminals of the battery, Once a short circuit is found ,it will cause the battery heat, fire and other dangers;

★Do not burn or heat the battery. Heating or burning the battery will melt the insulation, thus causing the battery fire, explosion and other dangers;

★Do not immerse the battery in water or in a humidity environment, The water in the battery pack is easy to cause the short circuit of the protective circuit board in the battery pack, Thus causing short circuit heating and other dangers;

★You must charge the battery with a charger that matches the battery, Otherwise, abnormal chemical reactions will occur in the battery, causing the battery to heat up, burn and other dangers;

★It is strictly prohibited to use a hammer or other hard objects to hit the battery, cutting the battery with sharp tools, When the battery is damaged, it may cause short circuit, thus causing fire

and other dangers;

★If the battery is damaged or leaked, and the leakage material enters the eyes carelessly, please don't rub it with your hands, wash it with clean water immediately, and go to the doctor as soon as possible;

★Do not mix lithium ion batteries made by different manufacturers or with different capacity. If so, it may cause heat, fire and other dangers in the process of charging and discharging;

★The battery needs to be fully charged in time (within 3 days) after use, If the battery pack is not used for a long time, it must be fully charged every two months, otherwise the battery will be damaged easily;

★Do not exceed the battery charging time (generally within 10 hours), for fear that the danger of overheating. If the battery charging exceeds the specified time and is not fully charged, please stop charging;

★Don't put the battery where a child can touch so as to avoid danger;

Do not store the battery in microwave or high pressure environment so as to avoid danger;

Attentions for safe driving

In order to ensure that you can use this portable electric folding car safely, please be sure to follow the following requirements for proper operation:

- ★Before driving for the first time, be sure to do driving exercises in the open field. Only drive on the road after mastering acceleration, deceleration, braking, honking and other performances perfectly
- ★Follow traffic rules and maintain moderate speed;Driving at top speed will not only dramatically increase power consumption and reduce range, but also reduce safety
- ★When driving in a crowded indoor environment, such as a supermarket, you must keep a low speed to control the braking distance.When turning, the speed must not exceed 2km/h
- ★It should not be used in rainy days, and wading should be avoided as far as possible. When wading through, the wading depth should not exceed 20mm and the speed should not exceed 1km/h
- ★With the power switch key, in the non-driving state, the power switch must be closed and the key pulled out, to prevent children and others from wrong operation

Eliminate of frequent malfunction

Malfunction	Reason	solution
Adjuster shifter out of order or Significant reduction of maximum speed	1. battery cable is loose 2. Shortage power 3. Spring failure of Adjuster shifter	1.Tighten the battery cable 2.full recharge the battery 3.change the spring or shifter
Charger cannot be charged normally	1.Power line damaged or socket off 2.Charger fuse fusing 3.The battery terminal is disconnected	1.Repair the power cord or tighten the socket 2.Replace fuse of the same model 3.Plug in all the connections
Charging is normal,but driving mileage is not enough	1. Tire pressure is insufficient or brake disc adjustment is too tight 2. Long time driving at high speed or travel uphill, adverse wind, brake start too frequently 3. Battery aging or damage	1. Properly inflate tires;Adjust the tightness of brake cable properly 2. Avoid long and high speed driving and frequent braking;It is normal for power consumption to increase in uphill and upwind conditions 3. Replace or repair batteries
drive motor doesn't work properly	1. Loose or detached connection of battery, motor, controller and other cables 2. Can not maintain the average speed of driving, control shifter is damage caused by pollution	1. Check the battery, motor, controller and other cables, and reconnect firmly 2.Keep the control shifter clean to avoid stained.

★Batteries, chargers, motors, controllers and control shifter shall be repaired by professional technicians or after-sales maintenance units

★The above data are for reference only. The technical data of this product series are subject to change without prior notice

★In actual use, due to different road conditions, climate, driving speed, charging time, mileage and other data are different

Friendly reminder

★When the scooter starts, please slow down the throttle to prevent the front wheel from skidding

★When the scooter climbing the slope , please lean forward to Keep the center of gravity steady and prevent scooter from turning over after climbing the slope

★When the scooter is turning, please slow down and avoid turning over

★in the normal using of the scooter, please check whether the screw loose or off periodically

★Please drive as little as possible when the road is slippery in rainy days

**Specification is based solely on actual products;
Subject to change without prior notice.**