


Transfer chair

User manual

Model: GE-HY803

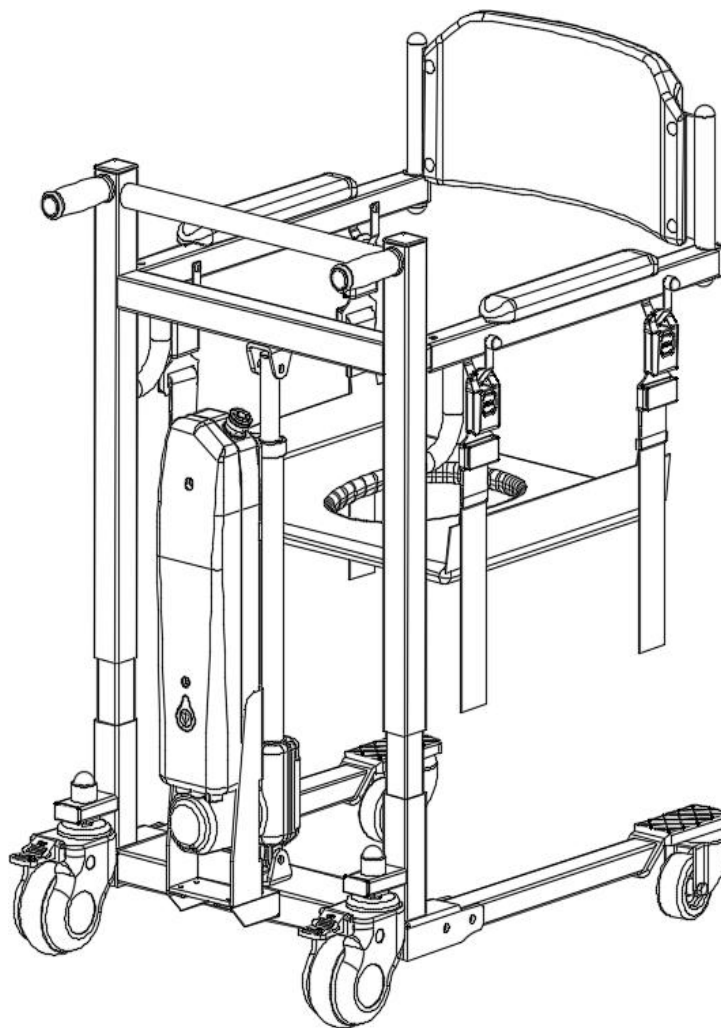


Gilani Engineering Pty Ltd

Contents

I. Parts introduction.....	1
II. Technical specification.....	1
III. Installation instructions.....	2
IV. User guide.....	4
V. Troubleshooting.....	6
VI. Warranty.....	6


I. Parts introduction

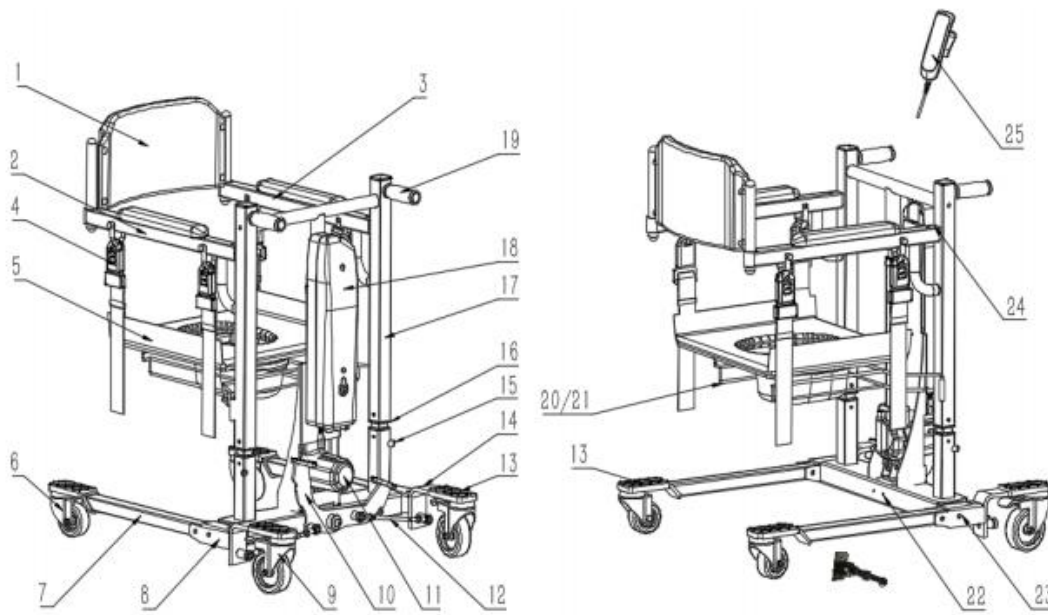

1、 Backrest; 2、 Right armrest; 3、 Left armrest; 4、 Ring; 5、 Cushion;
 6、 3" front wheel; 7、 Support leg; 8、 Right rear wheel support; 9、 4"Rear
 wheel; 10、 Pedal; 11、 Battery; 12、 Connecting rod; 13、 Castor mat;
 14、 Left rear wheel support; 15、 Screw; 16、 Sleeve; 17、 Lifting frame
 18、 Controller box/ Motor; 19、 Handle; 20/21; Bedpen/ Bedpen support; 22,
 Support; 23、 Long hinge; 24、 Short hinge; 25、 Controller


II. Technical specification

Overall size: 900x676x870mm	Rated power: 100~240Vac 50/60Hz 72VA
Net weight: 45 Kg	Voltage: DC 24V±1
Lifting range: 0~250mm	Storage temperature: —5~50℃
Lifting speed: 7.6mm/s	Work temperature: 5~40℃
Motor: DC24V actuator motor	Relative humidity: ≤85%RH
Charger: AC220V/50Hz	Barometric pressure: 86~106kPa
Loading capacity: 150Kg	Waterproof: IPX4
Noise: 65 (dB)	
Lifting times per full charge: 60 Times	

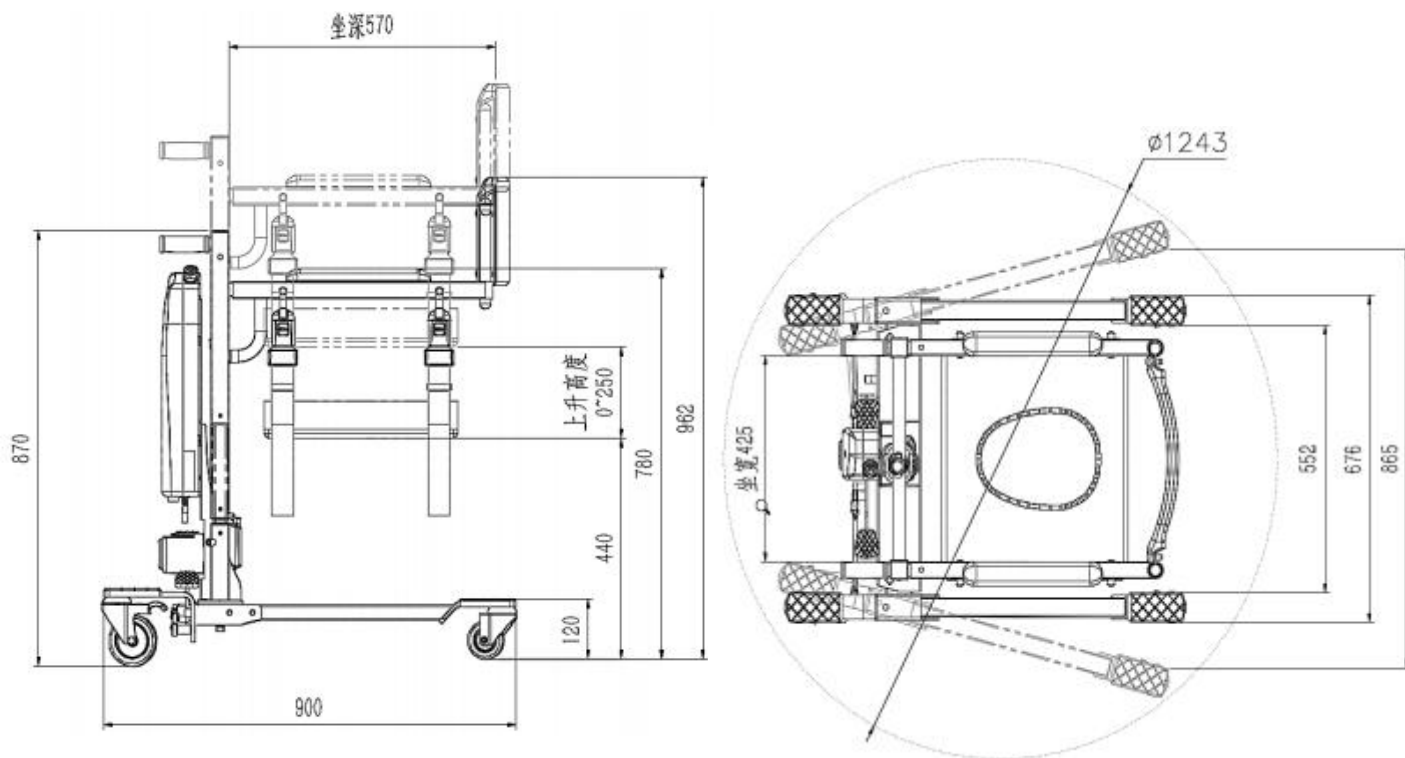

III. Installation instructions

1. Parts

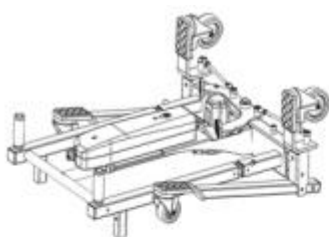

Frame x 1 set



Armrest 1 pair



Bedpan 1 pc




Charger cable



Long hinge x 2 pcs
Short hinge x 2 pcs




Backrest 1 pc

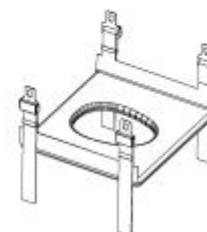


Cushion 1 pc

2. Installation steps


Step 1: install arm-brace

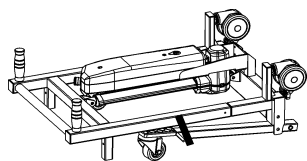


①Take out the body

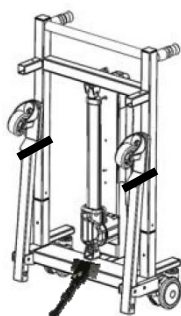


②unfold the arm-brace

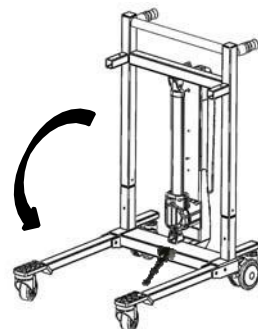


③Untie the strap, unfold the support leg

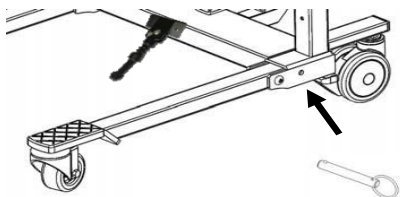

④Insert the long bolt into the hole

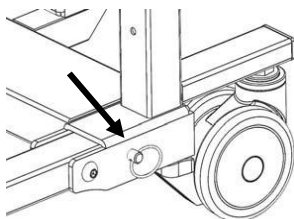


⑤Install bolts on both sides of the same operation

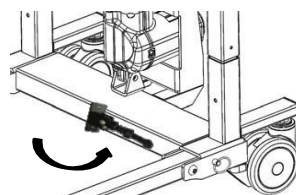


⑥Fold the arm-brace


Step 2: Install armrest

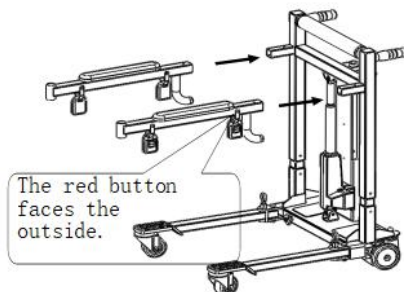


①Distinguish the left and right armrests.

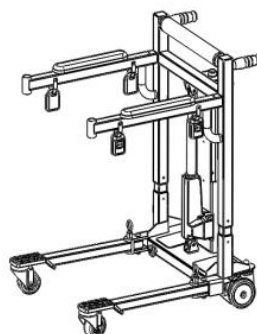

②Insert the two armrests into the main frame.

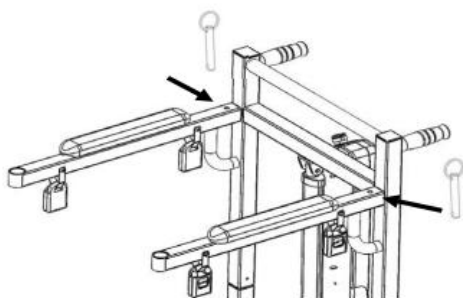

③ Insert the short bolt into the hole

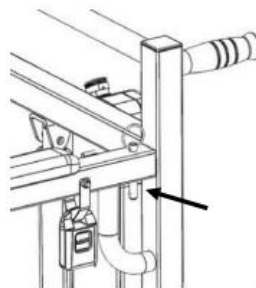


④ complet

IV. User guide

1. How to use the controller panel

Press the up and down button to control the machine, release the button to stop the movement.

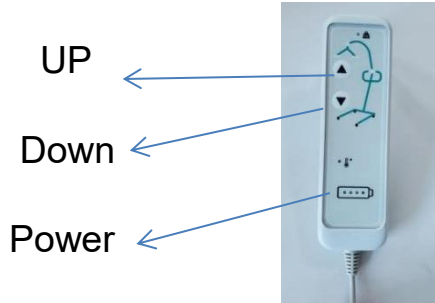



2. Emergency stop

The emergency stop button is used to interrupt and turn on the power. so that the motor will stop working immediately in the case of danger. Stop working immediately in case of danger.

(Close: Press the red button, stop working.①)


Open: Rotate the red button clockwise, normal working.②

Suggestion: If the machine has not been used for more than 3 days, push this button to protect the battery.

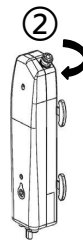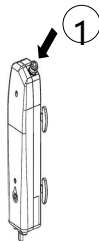

3. Control device connection





There were 4 sockets.

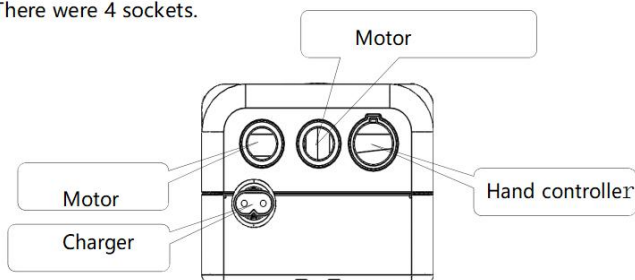

4. Front wheel brake




5. Operation of the transfer chair.

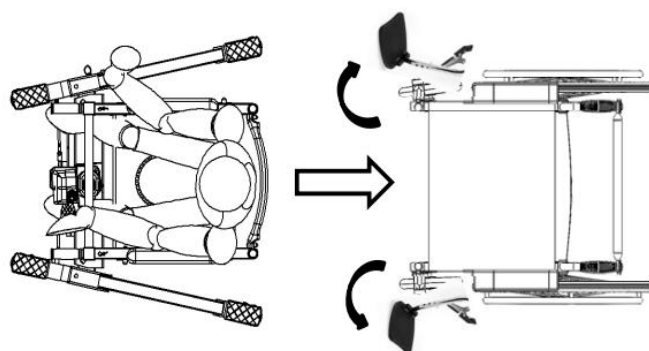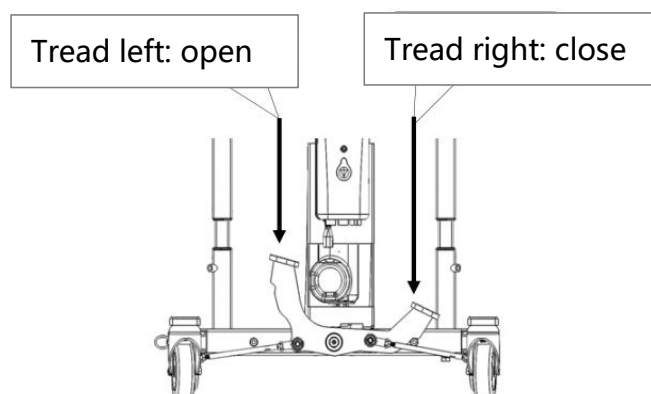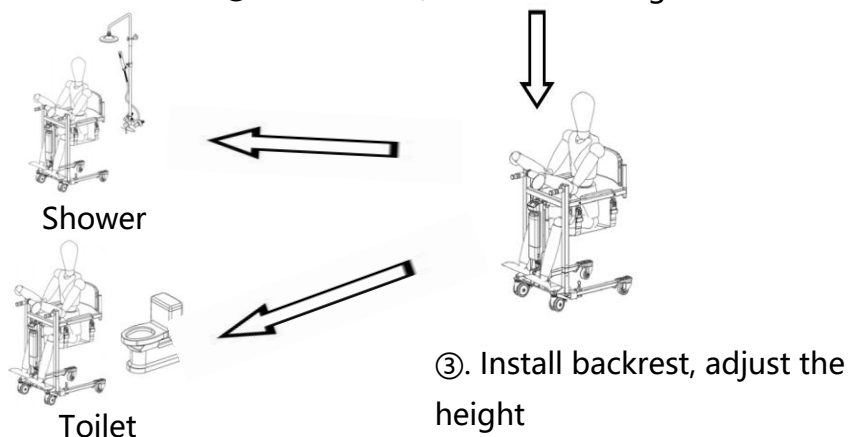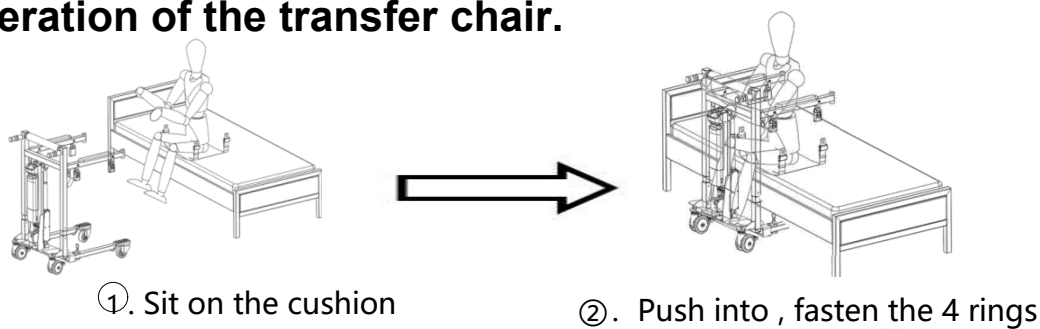

6. Storage of transfer chair

Folding steps:

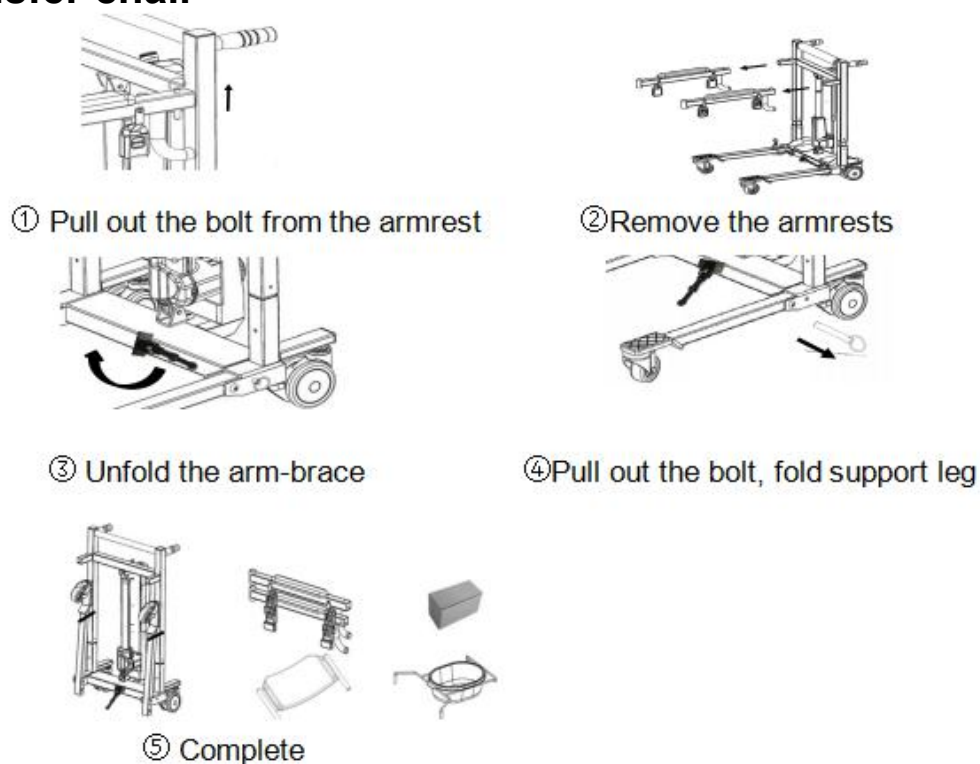

7. Battery charge.

When the battery power is only 1 grid, it needs to be charged immediately. If not, the motor will stop working.



Notes:

1. Please close the emergency button and turn off the machine when charging.
2. It is generally charged once every month when the machine is not used for a long time.
3. Do not charge in a humid environment.
4. If the shifter is not used every day, it is recommended to close the emergency stop button.

V. Troubleshooting

Couldn't lift or land	1. Check whether the emergency stop button is pressed 2. Check whether the wire of is connected correctly 3. Check battery power 4. Check whether the battery contact panel has problems or damages 5. If the electric shift machine is not working, please contact the after-sales service
Charging failure	1. Check whether the charger is plugged into the socket 2. Check the battery contact panel for problems or damage 3. If the charger still does not charge, please contact after-sales service
It stops while lifting	1. Check whether the emergency stop button is pressed 2. Use an emergency descent device 3. If the electric shift machine is not working, please contact the after-sales service.

VI. Warranty

Warranty Card			
User Name		Phone	
Address			
Product Name		Model	
Serial Number			
Maintenance Record	Date	Causes and Handing	Repairer