

MOBILITY ROBOT

M4/M4U/M4-Pro/M4U-Pro

Safety

- Avoid obstacles Automatically, Prevent to Fall
- No Skidding, No Hill-hopping
- Sufficient Power, Stable Speed
- Dual Mode Options

Comfort

- Automatic Postural Balance Control
- Wheelchair Height Adjustment

Convenience

- Automatic Fold/Unfold, Detachable
- Support APP Control



M4/M4-Pro

Mecanum Wheel Pedal

Smaller turn-in-place diameter for easier passage through narrow spaces



Maximum Load Capacity

115kg



Endurance Mileage

17km



Minimum Turning Radius

820mm



Weight

49.5kg



Maximum Speed

6km/h



Maximum Slope

15°



Maximum Load Capacity

115kg



Endurance Mileage

17km



Minimum Turning Radius

820mm



Weight

46.5kg



Maximum Speed

6km/h



Maximum Slope

8°

M4U/M4U-Pro

Omni Wheel Peda

Ultra-quiet, strong vibration damping, smooth operation



Powerful

The M4 is powered by a large-capacity battery and a high-performance motor, and is equipped with a closed-loop control system for smooth and safe operation. It can easily handle even rough roads and slopes.

Automatic Postural Balance Control

Based on the gyroscope balance technology, adding the self-developed advanced processing algorithm, M4 can adjust the tilt degree of the seat, keep the seat in a horizontal state to create a comfortable user experience





Folding Design, Height Adjustable

The M4 folds automatically and adjust to different heights. The drivers can fit different desks, such as office desks, dining tables, and bar counters. This makes it convenient for users to work and socialize.


Support APP Control

The M4 supports online control of mobile devices, APP operation is simpler and smarter. It is equipped with remote positioning, one-key summoning function.





MOBILITY ROBOT

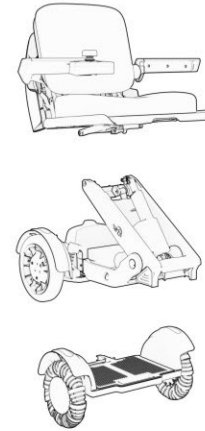

Modular Design

The M4 is composed of seat, executive unit and front wheels, which can be quickly assembled and disassembled. It is easy to carry, and can be stored in the trunk of the car.

Three-step quick Disassembly

The device can be easily disassembled and reassembled in three steps.



M4



M4-Pro



M4U



M4U-Pro



Function comparison

Product model	M4	M4U	M4-Pro	M4U-Pro
Self balancing seat position	✓	✓	✓	✓
Adjustable seat height	✓	✓	✓	✓
Folding and unfolding	✓	✓	✓	✓
Support APP control	✓	✓	✓	✓
Quick disassemble	✓	✓	✓	✓
Avoid obstacles automatically			✓	✓
Safety driving assist system			✓	✓
Intelligent protection in fall			✓	✓
Automatic driverless			✓	✓